

CERTIFICATE OF ANALYSIS FOR

ZINC CONCENTRATE

CERTIFIED REFERENCE MATERIAL

OREAS 350

Table 1. Certified Values, SDs, 95% Confidence and Tolerance Limits for OREAS 350.

Constituent	Certified Value	1SD	95% Confidence Limits		95% Tolerance Limits	
			Low	High	Low	High
Umpire Labs (dry sample basis)						
<i>Acid Digestion Titration*</i>						
Zn, Zinc (wt.%)	44.12	0.144	44.04	44.21	44.09	44.16
Commercial Labs ('as received' sample basis)						
<i>4-Acid Digestion</i>						
Ag, Silver (ppm)	96.6	2.68	95.5	97.7	94.6	98.6
As, Arsenic (ppm)	900	28	881	919	883	916
Bi, Bismuth (ppm)	< 5	IND	IND	IND	IND	IND
Cd, Cadmium (ppm)	932	41	909	955	917	947
Cu, Copper (wt.%)	0.854	0.014	0.847	0.862	0.839	0.869
Fe, Iron (wt.%)	5.23	0.182	5.12	5.33	5.15	5.31
Mg, Magnesium (wt.%)	0.263	0.018	0.253	0.273	0.258	0.269
Pb, Lead (wt.%)	8.22	0.211	8.08	8.35	8.09	8.34
Sb, Antimony (ppm)	108	20	93	123	104	112
Tl, Thallium (ppm)	128	11	121	134	125	130
Zn, Zinc (wt.%)	44.13	1.365	42.49	45.77	42.32	45.94
<i>Peroxide Fusion ICP</i>						
SiO ₂ , Silicon dioxide (wt.%)	6.12	0.232	5.99	6.26	5.99	6.26
<i>Infrared Combustion</i>						
S, Sulphur (wt.%)	28.82	1.093	28.17	29.48	28.36	29.28

*Except for one laboratory that used borate fusion XRF.
 Note: intervals may appear asymmetric due to rounding.

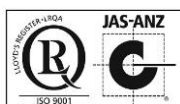


Table 2. Indicative Values for OREAS 350.

Constituent	Unit	Value	Constituent	Unit	Value	Constituent	Unit	Value
Umpire Labs (dry sample basis)								
<i>Acid Digestion Titration</i>								
Pb	wt.%	8.47						
<i>Thermogravimetry</i>								
H ₂ O-	wt.%	0.495						
Commercial Labs (as received basis)								
<i>4-Acid Digestion</i>								
Al	wt.%	0.640	La	ppm	5.00	Se	ppm	< 5
Ba	ppm	11.4	Li	ppm	6.00	Sn	ppm	< 10
Be	ppm	0.70	Mn	wt.%	0.023	Sr	ppm	4.00
Ca	wt.%	0.350	Mo	ppm	4.20	Te	ppm	< 5
Ce	ppm	< 10	Na	wt.%	0.030	Ti	wt.%	0.030
Co	ppm	24.0	Nb	ppm	< 10	U	ppm	< 10
Cr	ppm	21.0	Ni	ppm	12.0	V	ppm	15.4
Ga	ppm	13.0	P	wt.%	0.010	W	ppm	< 10
Ge	ppm	11.0	Re	ppm	< 5	Y	ppm	2.00
In	ppm	< 10	S	wt.%	28.86	Zr	ppm	13.0
K	wt.%	0.450	Sc	ppm	1.00			
<i>Peroxide Fusion ICP</i>								
Al	wt.%	0.596	Fe	wt.%	5.13	S	wt.%	28.63
As	ppm	835	K	wt.%	0.440	Sb	ppm	36.6
Ca	wt.%	0.380	Mg	wt.%	0.277	Sr	ppm	< 50
Cd	ppm	855	Mn	wt.%	0.023	Ti	wt.%	0.030
Co	ppm	25.0	Ni	ppm	< 50	Zn	wt.%	48.84
Cu	wt.%	0.846	Pb	wt.%	8.32			
<i>Thermogravimetry</i>								
H ₂ O-	wt.%	0.530						

Note: the number of significant figures reported is not a reflection of the level of certainty of stated values. They are instead an artefact of ORE's in-house CRM-specific LIMS.

INTRODUCTION

OREAS reference materials are intended to provide a low cost method of evaluating and improving the quality of analysis of geological samples. To the geologist they provide a means of implementing quality control in analytical data sets generated in exploration from the grass roots level through to prospect evaluation, and in grade control at mining operations. To the analyst they provide an effective means of calibrating analytical equipment, assessing new techniques and routinely monitoring in-house procedures.

SOURCE MATERIALS

OREAS 350 is a zinc concentrate certified reference material (CRM) prepared and certified by Ore Research & Exploration Pty Ltd. The material constituting OREAS 350 was sourced from the McArthur River Mining operation located in the Gulf of Carpentaria region in the Northern Territory of Australia. The mine is situated about 900 km southeast of Darwin. The major sulphides are pyrite, sphalerite and galena, with lesser chalcopyrite, arsenopyrite and marcasite.

COMMINUTION AND HOMOGENISATION PROCEDURES

The source materials constituting OREAS 350 were prepared in the following manner:

- drying of materials to constant mass at 85°C;
- deagglomeration and homogenisation;
- packaging into 10g units sealed under nitrogen in laminated foil pouches.

ANALYTICAL PROGRAM

Sixteen umpire laboratories and fourteen commercial analytical laboratories participated in the round robin certification program to certify the 14 analytes reported in Table 1.

The umpire laboratories each received a 30g sample and employed the following methods:

- Zinc was determined in 3 trials by ethylenediaminetetraacetic acid (EDTA) with titration finish (15 laboratories) or borate fusion with XRF finish (1 laboratory);
- Moisture (H₂O-) at 105°C by thermogravimetry (16 laboratories).

The umpire laboratories were given strict pre-assay sample instructions relating to moisture correction. These instructions included:

- Equilibration of sample material to lab atmosphere for a minimum of 2 hours;
- Hygroscopic moisture analysis at 105°C determined on a separate subsample and weighed for analysis at the same time as the sample aliquots for Zn as per ISO 9599.

The umpire laboratories were also requested to report metal concentrations on both a dry and moisture-bearing basis and include all results for moisture determinations. A value for moisture (H₂O-) is provided in the Indicative Values table for informational purposes only (see Table 2). **The certified value for zinc is on a dry sample basis.**

The commercial laboratories each received five 20g samples and employed the following methods:

- 4-acid digestion (HF-HNO₃-HClO₄-HCl) with ICP-OES, ICP-MS or AAS finish for Ag, As, Bi, Cd, Cu, Fe, Mg, Pb, Sb, Tl and Zn (14 laboratories);
- Peroxide fusion with ICP-OES finish for SiO₂ (14 laboratories);
- Infrared combustion furnace for S (13 laboratories).

Samples for the round robin program were taken at five predetermined sampling intervals immediately following final homogenisation and are considered representative of the entire batch of OREAS 350. The five samples received by each laboratory were obtained by taking one 20g sample from each of the five sampling lots. Table 1 presents the 14 certified values together with their associated 1SD's, 95% confidence and tolerance limits and Table 2 shows 52 indicative values. Table 3 provides performance gate intervals for the certified values based on their pooled 1SD's. Tabulated results of all elements together with uncorrected means, medians, standard deviations, relative standard deviations and per cent deviation of lab means from the corrected mean of means (PDM³) are presented in the detailed certification data for this CRM (**OREAS 350 DataPack.xlsx and OREAS 350 Umpire DataPack.xlsx**).

STATISTICAL ANALYSIS

Certified Values, Confidence Limits, Standard Deviations and Tolerance Limits (Table 1) have been determined for each analyte following removal of individual, laboratory dataset (batch) and 3SD outliers (single iteration). For individual outliers within a laboratory batch the z-score test is used in combination with a second method that determines the per cent deviation of the individual value from the batch median. Outliers in general are selected on the basis of z-scores > 2.5 and with per cent deviations (i) > 3 and (ii) more than three times the average absolute per cent deviation for the batch. In certain instances statistician's prerogative has been employed in discriminating outliers. Each laboratory data set mean is tested for outlying status based on z-score discrimination and rejected if >2.5 . After individual and laboratory data set (batch) outliers have been eliminated a non-iterative 3 standard deviation filter is applied, with those values lying outside this window also relegated to outlying status.

Certified Values are the means of accepted laboratory means after outlier filtering. Indicative (uncertified) values (Table 2) are provided where i) the number of laboratories reporting a particular analyte is insufficient (< 5) to support certification; ii) inter-laboratory consensus is poor; or iii) a significant proportion of results are outlying.

95% Confidence Limits are inversely proportional to the number of participating laboratories and inter-laboratory agreement. It is a measure of the reliability of the certified value. A 95% confidence interval indicates a 95% probability that the true value of the analyte under consideration lies between the upper and lower limits. *95% Confidence Limits should not be used as control limits for laboratory performance.*

Standard Deviation values (1SDs) are reported in Table 1 and provide an indication of a level of performance that might reasonably be expected from a laboratory being monitored by this CRM in a QA/QC program. The SD values include all sources of measurement uncertainty: between-lab variance, within-run variance (precision errors) and CRM variability. For an effective CRM the contribution of the latter should be negligible in comparison to measurement errors. OREAS reference materials have a level of homogeneity such that the observed variance from repeated analysis has its origin almost exclusively in the analytical process rather than the reference material itself.

The SD for each analyte's certified value is calculated from the same filtered data set used to determine the certified value, i.e. after removal of any individual, lab dataset (batch) and 3SD outliers (single iteration). These outliers can only be removed after the absolute homogeneity of the CRM has been independently established, i.e. the outliers must be confidently deemed to be analytical rather than arising from inhomogeneity of the CRM. **The standard deviation is then calculated for each analyte from the pooled accepted analyses generated from the certification program.**

In the application of SD's in monitoring performance it is important to note that not all laboratories function at the same level of proficiency and that different methods in use at a particular laboratory have differing levels of precision. Each laboratory has its own inherent SD (for a specific concentration level and analyte-method pair) based on the analytical process and this SD is not directly related to the round robin program.

The majority of data generated in the round robin program was produced by a selection of world class laboratories. The SD's thus generated are more constrained than those that would be produced across a randomly selected group of laboratories. To produce more

generally achievable SD's the 'pooled' SD's provided in this report include inter-lab bias. This 'one size fits all' approach may require revision at the discretion of the QC manager concerned following careful scrutiny of QC control charts.

Table 3 shows **Performance Gates** calculated for two and three standard deviations. As a guide these intervals may be regarded as warning or rejection for multiple 2SD outliers, or rejection for individual 3SD outliers in QC monitoring, although their precise application should be at the discretion of the QC manager concerned. A second method utilises a 5% window calculated directly from the certified value. Standard deviation is also shown in relative per cent for one, two and three relative standard deviations (1RSD, 2RSD and 3RSD) to facilitate an appreciation of the magnitude of these numbers and a comparison with the 5% window. Caution should be exercised when concentration levels approach lower limits of detection of the analytical methods employed as performance gates calculated from standard deviations tend to be excessively wide whereas those determined by the 5% method are too narrow.

Table 3. Performance Gates for OREAS 350.

Constituent	Certified Value	Absolute Standard Deviations					Relative Standard Deviations			5% window	
		1SD	2SD Low	2SD High	3SD Low	3SD High	1RSD	2RSD	3RSD	Low	High
Umpire Labs (dry sample basis)											
<i>Acid Digestion Titration*</i>											
Zn, wt.%	44.12	0.144	43.83	44.41	43.69	44.55	0.33%	0.65%	0.98%	41.92	46.33
Commercial Labs ('as received' sample basis)											
<i>4-Acid Digestion</i>											
Ag, ppm	96.6	2.68	91.3	102.0	88.6	104.6	2.77%	5.54%	8.31%	91.8	101.4
As, ppm	900	28	844	956	816	984	3.12%	6.24%	9.37%	855	945
Bi, ppm	< 5	IND	IND	IND	IND	IND	IND	IND	IND	IND	IND
Cd, ppm	932	41	849	1015	808	1056	4.44%	8.88%	13.32%	886	979
Cu, wt.%	0.854	0.014	0.826	0.882	0.812	0.896	1.65%	3.30%	4.95%	0.811	0.897
Fe, wt.%	5.23	0.182	4.86	5.59	4.68	5.78	3.49%	6.98%	10.47%	4.97	5.49
Mg, wt.%	0.263	0.018	0.228	0.299	0.210	0.317	6.79%	13.59%	20.38%	0.250	0.276
Pb, wt.%	8.22	0.211	7.79	8.64	7.58	8.85	2.57%	5.14%	7.70%	7.80	8.63
Sb, ppm	108	20	67	149	47	169	18.92%	37.83%	56.75%	102	113
Tl, ppm	128	11	106	149	96	160	8.36%	16.72%	25.08%	121	134
Zn, wt.%	44.13	1.365	41.40	46.86	40.03	48.23	3.09%	6.19%	9.28%	41.92	46.34
<i>Peroxide Fusion ICP</i>											
SiO ₂ , wt.%	6.12	0.232	5.66	6.59	5.43	6.82	3.79%	7.58%	11.36%	5.82	6.43
<i>Infrared Combustion</i>											
S, wt.%	28.82	1.093	26.64	31.01	25.54	32.10	3.79%	7.58%	11.38%	27.38	30.26

*Except for one laboratory that used borate fusion XRF.
Note: intervals may appear asymmetric due to rounding.

Tolerance Limits (ISO Guide 3207) were determined using an analysis of precision errors method and are considered a conservative estimate of true homogeneity. The meaning of tolerance limits may be illustrated for Zn by the umpire laboratories employing EDTA digestion with titration, where 99% of the time ($1-\alpha=0.99$) at least 95% of subsamples ($\rho=0.95$) will have concentrations lying between 44.09 and 44.16 wt.%. Put more precisely, this means that if the same number of subsamples were taken and analysed in the same manner repeatedly, 99% of the tolerance intervals so constructed would cover at least 95% of the total population, and 1% of the tolerance intervals would cover less than 95% of the total population (ISO Guide 35). Based on the statistical analysis of the results of the inter-laboratory certification program it can be concluded that OREAS 350 is fit-for-purpose as a certified reference material (see 'Intended Use' below).

PARTICIPATING LABORATORIES

Actlabs, Ancaster, Ontario, Canada
AH Knight, St Helens, Merseyside, UK
Alex Stewart International, Liverpool, UK
ALS Inspection, Prescott, Merseyside, UK
ALS, Burnie, TAS, Australia
ALS, Lima, Peru
ALS, Perth, WA, Australia
ALS, Vancouver, BC, Canada
Bachelet, Angleur, Liege, Belgium
Bureau Veritas Geoanalytical, Adelaide, SA, Australia
Bureau Veritas Geoanalytical, Perth, WA, Australia
Independent, Perth, WA, Australia
Inspectorate (BV), Lima, Peru
Inspectorate (BV), Shanghai, Bao Shan District, China
Inspectorate (BV), Witham, Essex, UK
Intertek Genalysis, Perth, WA, Australia
Intertek LSI, Rotterdam, Zuid-Holland, Netherlands
Intertek Testing Services, Cupang, Muntinlupa, Philippines
MinAnalytical Services, Perth, WA, Australia
RC Inspection, Rotterdam, Netherlands
SGS Australia Mineral Services, Perth, WA, Australia
SGS Mineral Services, Townsville, QLD, Australia
SGS Nederland B.V., Spijkenisse, Zuid-Holland, Netherlands
Shiva Analyticals Ltd, Bangalore North, Karnataka, India
SRL, Perth, WA, Australia
Actlabs, Ancaster, Ontario, Canada
AH Knight, St Helens, Merseyside, UK
Alex Stewart International, Liverpool, UK

PREPARER AND SUPPLIER OF THE REFERENCE MATERIAL

Reference material OREAS 350 has been prepared, certified and is supplied by:

ORE Research & Exploration Pty Ltd
37A Hosie Street
Bayswater North VIC 3153
AUSTRALIA

Tel: +613-9729 0333
Fax: +613-9729 8338
Web: www.ore.com.au
Email: info@ore.com.au

It is available in unit sizes of 10g sealed under nitrogen in laminated foil pouches.

INTENDED USE

OREAS 350 is intended for the following uses:

- for the monitoring of laboratory performance in the analysis of analytes reported in Table 1 in geological samples;
- for the verification of analytical methods for analytes reported in Table 1;
- for the calibration of instruments used in the determination of the concentration of analytes reported in Table 1.

STABILITY AND STORAGE INSTRUCTIONS

OREAS 350 is a sulphide-rich reference material (S ~28.8%) and is reactive under normal atmospheric conditions. To inhibit oxidation and prolong its shelf life it has been sealed under nitrogen in robust laminated foil pouches. In its unopened state under normal conditions of storage it has a shelf life beyond five years.

INSTRUCTIONS FOR CORRECT USE

The certified value for Zn by EDTA with titration is on a dry sample basis. There was considerable variation in moisture content reported by the laboratories (lab means varied from 0.21 to 0.99 wt.% H₂O-) and this can be a significant source of error if not properly controlled. Therefore, all analyses were performed on the samples after equilibration with the laboratory atmosphere for a minimum of 2 hours and hygroscopic moisture analysis at 105°C determined on a separate subsample and weighed for analysis at the same time as the sample aliquot for Zn as per ISO 9599. If the reference material is not dried prior to analysis, the EDTA with titration certified value for Zn should be corrected to the moisture-bearing basis.

As per routine analysis at commercial laboratories, the certified values derived by 4-acid digestion, peroxide fusion with ICP-OES (SiO₂ only) and infrared combustion furnace (S only) refer to the concentration levels in the packaged state. There is no need for drying prior to weighing and analysis.

TRACEABILITY

The analytical samples were selected in a manner to represent the entire batch of prepared CRM. This 'representivity' was maintained in each submitted laboratory sample batch and ensures the user that the data is traceable from sample selection through to the analytical results that underlie the consensus values. Each analytical data set has been validated by its assayer through the inclusion of internal reference materials and QC checks during analysis. The laboratories were chosen on the basis of their competence (from past performance in inter-laboratory programs) for a particular analytical method, analyte, or analyte suite, and sample matrix. Most of these laboratories have and maintain ISO 17025 accreditation. The certified and non-certified (indicative) values presented in this report are calculated from the means of accepted data following robust statistical treatment as detailed in this report.

HANDLING INSTRUCTIONS

Fine powders pose a risk to eyes and lungs and therefore standard precautions such as the use of safety glasses and dust masks are advised.

LEGAL NOTICE

Ore Research & Exploration Pty Ltd has prepared and statistically evaluated the property values of this reference material to the best of its ability. The Purchaser by receipt hereof releases and indemnifies Ore Research & Exploration Pty Ltd from and against all liability and costs arising from the use of this material and information.

QMS ACCREDITED

ORE Pty Ltd is accredited to ISO 9001:2008 by Lloyd's Register Quality Assurance Ltd for its quality management system including development, manufacturing, certification and supply of CRMs.



CERTIFYING OFFICER

A handwritten signature in blue ink, appearing to read 'Craig Hamlyn', is positioned above a horizontal line.

Craig Hamlyn (B.Sc. Hons - Geology), Technical Manager - ORE P/L

REFERENCES

ISO Guide 30 (1992), Terms and definitions used in connection with reference materials.

ISO Guide 31 (2000), Reference materials – Contents of certificates and labels.

ISO Guide 3207 (1975), Statistical interpretation of data - Determination of a statistical tolerance interval.

ISO Guide 35 (2006), Certification of reference materials - General and statistical principals.