

CERTIFICATE OF ANALYSIS FOR
CERTIFIED REFERENCE MATERIAL
OREAS 250c
Gold Oxide Ores (Orogenic Gold Lodes, Australia)

Table 1. Certified Value, Uncertainty & Tolerance Intervals for Au by FA in OREAS 250c.

Constituent	Certified Value [†]	95% Expanded Uncertainty		95% Tolerance Limits	
		Low	High	Low	High
Pb Fire Assay					
Au, Gold (ppm)	0.313	0.306	0.319	0.310*	0.315*

SI unit equivalents: ppm (parts per million; 1×10^{-6}) \equiv mg/kg.


[†]This operationally defined measurand meets the requirements of ISO 17034 and all participating laboratories comply with the requirements of ISO 17025.

*Gold Tolerance Limits for typical 30g fire assay are determined from 20 x 85mg INAA results and the Sampling Constant (Ingamells & Switzer, 1973).

Note: intervals may appear asymmetric due to rounding.



Accredited for compliance with ISO 17034

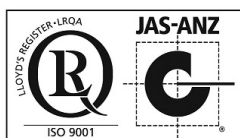

COA-1680-OREAS250c-R2 BUP-70-10-01 Ver:2.0	18-Mar-2024
---	-------------

Table 2. Certified Value, Uncertainty & Tolerance Intervals for other measurands in OREAS 250c.

Constituent	Certified Value	95% Expanded Uncertainty		95% Tolerance Limits	
		Low	High	Low	High
Aqua Regia Digestion (sample weights 10-50g)					
Au, Gold (ppm)	0.313	0.305	0.321	0.310*	0.316*
Cyanide Leach					
Au, Gold (ppm)	0.295	0.287	0.303	0.294*	0.296*
PhotonAssay					
Au, Gold (ppm)	0.311	0.302	0.320	0.310*	0.312*
4-Acid Digestion					
Ag, Silver (ppm)	< 0.1	IND	IND	IND	IND
Al, Aluminium (wt.%)	5.18	5.04	5.31	5.10	5.25
As, Arsenic (ppm)	54	51	57	52	56
Ba, Barium (ppm)	354	342	366	346	363
Be, Beryllium (ppm)	1.97	1.89	2.05	1.89	2.05
Bi, Bismuth (ppm)	0.29	0.26	0.32	0.27	0.31
Ca, Calcium (wt.%)	1.72	1.66	1.77	1.68	1.75
Cd, Cadmium (ppm)	0.075	0.044	0.105	0.064	0.085
Ce, Cerium (ppm)	70	67	74	67	74
Co, Cobalt (ppm)	20.6	19.5	21.7	20.0	21.1
Cr, Chromium (ppm)	100	95	106	96	104
Cs, Caesium (ppm)	2.62	2.48	2.75	2.52	2.71
Cu, Copper (ppm)	24.1	22.7	25.6	23.2	25.1
Dy, Dysprosium (ppm)	3.67	3.34	3.99	3.44	3.89
Er, Erbium (ppm)	1.59	1.44	1.75	1.50	1.69
Eu, Europium (ppm)	1.60	1.45	1.75	1.51	1.69
Fe, Iron (wt.%)	4.04	3.90	4.17	3.94	4.13
Ga, Gallium (ppm)	15.3	14.2	16.4	14.8	15.8
Gd, Gadolinium (ppm)	5.20	4.67	5.74	5.01	5.40
Hf, Hafnium (ppm)	4.34	4.08	4.61	4.11	4.58
Ho, Holmium (ppm)	0.63	0.57	0.68	0.59	0.67
In, Indium (ppm)	0.059	0.049	0.068	0.053	0.064
K, Potassium (wt.%)	1.32	1.27	1.36	1.29	1.35
La, Lanthanum (ppm)	38.0	35.8	40.2	36.8	39.2
Li, Lithium (ppm)	21.7	20.7	22.8	21.0	22.5
Lu, Lutetium (ppm)	0.18	0.14	0.21	IND	IND
Mg, Magnesium (wt.%)	1.55	1.50	1.61	1.52	1.59
Mn, Manganese (wt.%)	0.049	0.047	0.051	0.048	0.050
Mo, Molybdenum (ppm)	1.96	1.81	2.12	1.90	2.03
Na, Sodium (wt.%)	1.07	1.04	1.11	1.05	1.10

SI unit equivalents: ppm (parts per million; 1×10^{-6}) \equiv mg/kg; wt.% (weight per cent) \equiv % (mass fraction).

*Gold Tolerance Limits for typical 25g aqua regia digestion and 200g cyanide leach methods and 350g PhotonAssay methods are determined from 20 x 85mg INAA results and the Sampling Constant (Ingamells & Switzer, 1973).

Note: intervals may appear asymmetric due to rounding; IND = indeterminate (due to limited reading resolution of the methods employed. For practical purposes the 95% Expanded Uncertainty can be set between zero and a two times multiple of the upper bound/non-detect limit value).

Table 2 continued.

Constituent	Certified Value	95% Expanded Uncertainty		95% Tolerance Limits	
		Low	High	Low	High
4-Acid Digestion continued					
Nb, Niobium (ppm)	31.8	30.1	33.5	30.6	33.0
Nd, Neodymium (ppm)	32.1	29.0	35.2	31.2	33.0
Ni, Nickel (ppm)	83	80	87	81	85
P, Phosphorus (wt.%)	0.095	0.092	0.099	0.093	0.097
Pb, Lead (ppm)	9.92	9.09	10.75	9.41	10.43
Pr, Praseodymium (ppm)	8.51	7.73	9.29	8.10	8.93
Rb, Rubidium (ppm)	65	61	68	63	66
Re, Rhenium (ppm)	< 0.002	IND	IND	IND	IND
S, Sulphur (wt.%)	0.027	0.025	0.030	0.026	0.028
Sb, Antimony (ppm)	0.88	0.81	0.96	0.84	0.92
Sc, Scandium (ppm)	9.89	9.33	10.45	9.63	10.15
Sm, Samarium (ppm)	6.27	5.82	6.72	5.97	6.57
Sn, Tin (ppm)	3.59	3.37	3.82	3.34	3.84
Sr, Strontium (ppm)	311	300	323	304	319
Ta, Tantalum (ppm)	2.07	1.92	2.23	1.95	2.20
Tb, Terbium (ppm)	0.69	0.65	0.73	0.67	0.70
Th, Thorium (ppm)	8.68	8.23	9.12	8.34	9.01
Ti, Titanium (wt.%)	0.567	0.548	0.586	0.557	0.577
Tl, Thallium (ppm)	0.35	0.33	0.38	0.33	0.37
Tm, Thulium (ppm)	0.20	0.18	0.22	0.19	0.21
U, Uranium (ppm)	1.62	1.51	1.73	1.55	1.69
V, Vanadium (ppm)	70	67	72	68	72
W, Tungsten (ppm)	2.05	1.90	2.21	1.90	2.21
Y, Yttrium (ppm)	15.7	14.6	16.7	15.1	16.2
Yb, Ytterbium (ppm)	1.28	1.15	1.42	1.19	1.37
Zn, Zinc (ppm)	75	72	78	73	77
Zr, Zirconium (ppm)	188	179	196	184	192
Aqua Regia Digestion					
Ag, Silver (ppm)	0.079	0.067	0.091	IND	IND
Al, Aluminium (wt.%)	1.64	1.56	1.73	1.59	1.69
As, Arsenic (ppm)	43.4	42.0	44.9	42.6	44.3
B, Boron (ppm)	< 10	IND	IND	IND	IND
Ba, Barium (ppm)	95	89	100	91	98
Be, Beryllium (ppm)	1.08	1.01	1.14	1.04	1.11
Bi, Bismuth (ppm)	0.26	0.25	0.28	0.25	0.28
Ca, Calcium (wt.%)	0.524	0.493	0.555	0.511	0.538
Cd, Cadmium (ppm)	0.067	0.048	0.087	IND	IND

SI unit equivalents: ppm (parts per million; 1×10^{-6}) \equiv mg/kg; wt.% (weight per cent) \equiv % (mass fraction).

Note: intervals may appear asymmetric due to rounding; IND = indeterminate (due to limited reading resolution of the methods employed. For practical purposes the 95% Expanded Uncertainty can be set between zero and a two times multiple of the upper bound/non-detect limit value).

Table 2 continued.

Constituent	Certified Value	95% Expanded Uncertainty		95% Tolerance Limits	
		Low	High	Low	High
Aqua Regia Digestion continued					
Ce, Cerium (ppm)	40.9	39.2	42.5	39.7	42.0
Co, Cobalt (ppm)	15.9	15.1	16.7	15.1	16.8
Cr, Chromium (ppm)	45.9	42.4	49.4	44.5	47.3
Cs, Caesium (ppm)	1.13	1.05	1.20	1.08	1.17
Cu, Copper (ppm)	18.9	17.6	20.2	18.1	19.7
Fe, Iron (wt.%)	3.01	2.85	3.16	2.95	3.07
Ga, Gallium (ppm)	5.97	5.49	6.44	5.78	6.15
Ge, Germanium (ppm)	0.100	0.076	0.123	IND	IND
In, Indium (ppm)	0.025	0.019	0.032	0.022	0.029
K, Potassium (wt.%)	0.335	0.317	0.353	0.325	0.345
La, Lanthanum (ppm)	23.2	22.3	24.0	22.5	23.8
Li, Lithium (ppm)	7.61	7.12	8.10	7.32	7.90
Mg, Magnesium (wt.%)	0.989	0.942	1.035	0.967	1.010
Mn, Manganese (wt.%)	0.033	0.032	0.035	0.033	0.034
Mo, Molybdenum (ppm)	1.32	1.23	1.42	1.28	1.37
Na, Sodium (wt.%)	0.351	0.331	0.370	0.338	0.363
Nb, Niobium (ppm)	0.59	0.48	0.69	0.53	0.64
Ni, Nickel (ppm)	73	70	75	71	75
P, Phosphorus (wt.%)	0.058	0.056	0.061	0.057	0.060
Pb, Lead (ppm)	7.08	6.66	7.50	6.81	7.36
Rb, Rubidium (ppm)	19.1	18.1	20.1	18.5	19.7
Re, Rhenium (ppm)	< 0.001	IND	IND	IND	IND
S, Sulphur (wt.%)	0.024	0.020	0.028	0.023	0.026
Sb, Antimony (ppm)	0.45	0.40	0.49	0.42	0.47
Sc, Scandium (ppm)	3.74	3.51	3.97	3.60	3.87
Sn, Tin (ppm)	1.24	1.12	1.36	IND	IND
Sr, Strontium (ppm)	63	61	66	62	65
Ta, Tantalum (ppm)	< 0.01	IND	IND	IND	IND
Tb, Terbium (ppm)	0.45	0.40	0.49	IND	IND
Th, Thorium (ppm)	5.85	5.49	6.21	5.66	6.04
Ti, Titanium (wt.%)	0.162	0.151	0.174	0.157	0.167
Tl, Thallium (ppm)	0.14	0.12	0.15	IND	IND
U, Uranium (ppm)	0.97	0.89	1.04	0.94	0.99
V, Vanadium (ppm)	33.7	32.1	35.3	32.5	34.9
W, Tungsten (ppm)	0.38	0.33	0.44	0.33	0.44
Y, Yttrium (ppm)	9.10	8.74	9.46	8.90	9.30
Yb, Ytterbium (ppm)	0.63	0.55	0.71	IND	IND

SI unit equivalents: ppm (parts per million; 1×10^{-6}) \equiv mg/kg; wt.% (weight per cent) \equiv % (mass fraction).

Note: intervals may appear asymmetric due to rounding; IND = indeterminate (due to limited reading resolution of the methods employed. For practical purposes the 95% Expanded Uncertainty can be set between zero and a two times multiple of the upper bound/non-detect limit value).

Table 2 continued.

Constituent	Certified Value	95% Expanded Uncertainty		95% Tolerance Limits	
		Low	High	Low	High
Aqua Regia Digestion continued					
Zn, Zinc (ppm)	46.4	43.8	49.0	45.3	47.5
Zr, Zirconium (ppm)	37.1	28.7	45.6	35.2	39.1

SI unit equivalents: ppm (parts per million; 1×10^{-6}) \equiv mg/kg; wt.% (weight per cent) \equiv % (mass fraction).

Note: intervals may appear asymmetric due to rounding; IND = indeterminate (due to limited reading resolution of the methods employed).

Table 3. Indicative Values for OREAS 250c.

Constituent	Unit	Value	Constituent	Unit	Value	Constituent	Unit	Value
4-Acid Digestion								
Ge	ppm	0.14	Se	ppm	0.76			
Hg	ppm	0.030	Te	ppm	0.029			
Aqua Regia Digestion								
Dy	ppm	2.24	Ho	ppm	0.36	Se	ppm	0.19
Er	ppm	0.87	Lu	ppm	0.090	Sm	ppm	4.00
Eu	ppm	0.97	Nd	ppm	22.2	Te	ppm	0.023
Gd	ppm	3.38	Pd	ppb	< 10	Tm	ppm	0.11
Hf	ppm	0.68	Pr	ppm	5.77			
Hg	ppm	0.022	Pt	ppb	< 5			
Borate Fusion XRF								
Al ₂ O ₃	wt.%	10.07	MgO	wt.%	2.67	SiO ₂	wt.%	72.25
CaO	wt.%	2.35	MnO	wt.%	0.064	SO ₃	wt.%	0.070
Fe ₂ O ₃	wt.%	5.90	Na ₂ O	wt.%	1.46	TiO ₂	wt.%	0.997
K ₂ O	wt.%	1.58	P ₂ O ₅	wt.%	0.219			
Thermogravimetry								
LOI ¹⁰⁰⁰	wt.%	2.16						
Infrared Combustion								
C	wt.%	0.075	S	wt.%	0.010			
Laser Ablation ICP-MS								
Ag	ppm	0.100	Hf	ppm	5.99	Sm	ppm	6.32
As	ppm	51	Ho	ppm	0.77	Sn	ppm	4.50
Ba	ppm	350	In	ppm	0.038	Sr	ppm	296
Be	ppm	2.00	La	ppm	38.3	Ta	ppm	2.21
Bi	ppm	0.30	Lu	ppm	0.23	Tb	ppm	0.77
Cd	ppm	< 0.1	Mn	wt.%	0.051	Te	ppm	< 0.2
Ce	ppm	68	Mo	ppm	1.80	Th	ppm	9.09
Co	ppm	20.5	Nb	ppm	32.6	Ti	wt.%	0.599
Cr	ppm	115	Nd	ppm	33.1	Tl	ppm	0.15
Cs	ppm	2.52	Ni	ppm	86	Tm	ppm	0.29
Cu	ppm	32.0	Pb	ppm	10.5	U	ppm	1.81
Dy	ppm	4.12	Pr	ppm	8.90	V	ppm	73

SI unit equivalents: ppb (parts per billion; 1×10^{-9}) \equiv μ g/kg; ppm (parts per million; 1×10^{-6}) \equiv mg/kg; wt.% (weight per cent) \equiv % (mass fraction).

Note: the number of significant figures reported is not a reflection of the level of certainty of stated values. They are instead an artefact of ORE's in-house CRM-specific LIMS.

Table 3 continued.

Constituent	Unit	Value	Constituent	Unit	Value	Constituent	Unit	Value
Laser Ablation ICP-MS continued								
Er	ppm	2.04	Rb	ppm	63	W	ppm	2.50
Eu	ppm	1.59	Re	ppm	< 0.01	Y	ppm	19.8
Ga	ppm	14.6	Sb	ppm	0.90	Yb	ppm	1.75
Gd	ppm	5.33	Sc	ppm	9.55	Zn	ppm	75
Ge	ppm	1.08	Se	ppm	< 5	Zr	ppm	243

SI unit equivalents: ppm (parts per million; 1×10^{-6}) \equiv mg/kg.

Note: the number of significant figures reported is not a reflection of the level of certainty of stated values. They are instead an artefact of ORE's in-house CRM-specific LIMS.

TABLE OF CONTENTS

INTRODUCTION	8
SOURCE MATERIAL	8
COMMUNITION AND HOMOGENISATION PROCEDURES	9
PHYSICAL PROPERTIES.....	9
MINERALOGY	9
ANALYTICAL PROGRAM.....	10
STATISTICAL ANALYSIS	10
Homogeneity Evaluation.....	11
PERFORMANCE GATES	13
PARTICIPATING LABORATORIES.....	17
PREPARER AND SUPPLIER	17
METROLOGICAL TRACEABILITY	22
COMMUTABILITY	22
INTENDED USE.....	23
MINIMUM SAMPLE SIZE.....	23
PERIOD OF VALIDITY & STORAGE INSTRUCTIONS	23
INSTRUCTIONS FOR HANDLING & CORRECT USE	24
LEGAL NOTICE.....	25
DOCUMENT HISTORY.....	25
CERTIFYING OFFICER.....	25
QMS CERTIFICATION.....	25
REFERENCES	26

LIST OF TABLES

Table 1. Certified Value, Uncertainty & Tolerance Intervals for Au by FA in OREAS 250c.....	1
Table 2. Certified Value, Uncertainty & Tolerance Intervals for other measurands in OREAS 250c.....	2
Table 3. Indicative Values for OREAS 250c.....	5
Table 4. Physical properties of OREAS 250c.....	9
Table 5. Indicative mineralogy of OREAS 250c based on semi-quantitative XRD analysis.....	9
Table 6. Neutron Activation Analysis of Au on 20 x 85mg subsamples.....	13
Table 7. Performance Gates for OREAS 250c.....	14

LIST OF FIGURES

Figure 1. Au by Fire Assay in OREAS 250c	18
Figure 2. Au by Aqua Regia digestion in OREAS 250c.....	19
Figure 3. Au by Cyanide Leach in OREAS 250c	20
Figure 4. Au by PhotonAssay in OREAS 250c	21

INTRODUCTION

OREAS reference materials are intended to provide a low-cost method of evaluating and improving the quality of analysis of geological samples. To the geologist they provide a means of implementing quality control in analytical data sets generated in exploration from the grass roots level through to prospect evaluation, and in grade control at mining operations. To the analyst they provide an effective means of calibrating analytical equipment, assessing new techniques and routinely monitoring in-house procedures. OREAS reference materials enable users to successfully achieve process control of these tasks because the observed variance from repeated analysis has its origin almost exclusively in the analytical process rather than the reference material itself. In evaluating laboratory performance with this CRM, the section headed 'Instructions for correct use' should be read carefully.

Table 1 (all laboratories accredited to ISO 17025) and Table 2 (most laboratories accredited to ISO 17025) provide the certified values and their associated 95% expanded uncertainty and tolerance intervals, Table 3 shows indicative values including major and trace element characterisation, Table 4 provides some indicative physical properties and Table 5 provides indicative mineralogy based on semi-quantitative XRD analysis. Gold homogeneity (via INAA) is shown in Table 6 and is also demonstrated by a nested ANOVA (see 'Homogeneity Evaluation' section) and Table 7 presents the performance gate intervals for all certified values.

Tabulated results of all analytes together with uncorrected means, medians, standard deviations, relative standard deviations and per cent deviation of lab means from the corrected mean of means (PDM³) are presented in the detailed certification data for this CRM (**OREAS 250c-DataPack.1.2.240304_175341.xlsx**).

Results are also presented in scatter plots for gold by fire assay, aqua regia digestion, cyanide leach and PhotonAssay (Figures 1 to 4, respectively) together with $\pm 3SD$ (magenta) and $\pm 5\%$ (yellow) control lines and certified value (green line). Accepted individual results are coloured blue and individual and dataset outliers are identified in red and violet, respectively.

SOURCE MATERIAL

OREAS 250c was prepared from a blend of gold-bearing oxide ores and barren materials (basaltic scoria, quartz and mudstone). The ores were sourced from both the Magdala ore deposit (Stawell Gold Mine) located in west-central Victoria, Australia and the Wilber Lode (Andy Well Mine) located approximately 45km north of Meekatharra in the Murchison region of Western Australia.

Primary mineralisation in the Magdala ore consists of a quartz-sericite-carbonate schist assemblage containing the sulphides arsenopyrite, pyrrhotite and pyrite.

The Wilber Lode contains a common primary mineral assemblage of quartz, calcite, chlorite, fuchsite, pyrite, galena, sphalerite, chalcocopyrite and gold.

COMMINUTION AND HOMOGENISATION PROCEDURES

The material constituting OREAS 250c was prepared in the following manner:

- Drying of ore and barren materials to constant mass at 105°C;
- Crushing and milling of the barren materials to >98% minus 75 microns;
- Crushing and milling of the ore material to 100% minus 30 microns;
- Check analysis of ores for contained gold concentration;
- Blending ores and barren materials in appropriate proportions to achieve the desired grade;
- Homogenisation using OREAS' novel processing technologies;
- Packaging in 60g units sealed in laminated foil pouches and 500g units in plastic jars.

PHYSICAL PROPERTIES

OREAS 250c was tested at ORE Research & Exploration Pty Ltd's onsite facility for various physical properties. Table 4 presents these findings that should be used for informational purposes only.

Table 4. Physical properties of OREAS 250c.

Bulk Density (kg/m ³)	Moisture (wt.%)	Munsell Notation [‡]	Munsell Color [‡]
760	1.06	5YR 6/4	Light Brown

[‡]The Munsell Rock Color Chart helps geologists and archeologists communicate with colour more effectively by cross-referencing ISCC-NBS colour names with unique Munsell alpha-numeric colour notations for rock colour samples.

MINERALOGY

The semi-quantitative XRD results shown in Table 5 below have been normalised to 100 % and represent the relative proportion of crystalline material. Totals greater or less than 100 % are due to rounding errors. A trace amount of orthopyroxene might be present. 'Clay mineral' appears to be mainly smectite and/or vermiculite. Some amorphous material might be present.

Table 5. Indicative mineralogy of OREAS 250c based on semi-quantitative XRD analysis.

Mineral / Mineral Group	% (mass ratio)
Hematite	1
Goethite	0
Zeolite	0
Clay mineral	< 1
Kaolinite	2
Chlorite	< 1
Annite - biotite - phlogopite	1
Muscovite	6
Calcic amphibole	0
Clinopyroxene	4
Olivine	1
Tourmaline	0
Plagioclase	7
K-feldspar and/or rutile	1
Quartz	76
Dolomite - ankerite	0
Calcite	0

ANALYTICAL PROGRAM

Twenty-one commercial analytical laboratories participated in the program to certify the elements reported in Tables 1 and 2. The following methods were employed:

- Gold by fire assay (25-50g charge weight) with AAS (12 laboratories), ICP-OES (7 laboratories) finish and ICP-MS (2 laboratories) finish;
- Gold by aqua regia digestion (10-50g sample weight) with ICP-OES and/or ICP-MS (10 laboratories) finish and AAS finish (6 laboratories);
- Gold by cyanide leach; a variety of cyanide leach methods were undertaken by the participating laboratories including the use of LeachWELL tablets, alkaline added sodium cyanide solution as well as sodium cyanide liquor with LeachWELL powder. The sample weights included: 5g (2 laboratories by AAS finish), 20g (1 laboratory by AAS finish), 25g (1 laboratory by ICP-MS finish), 30g (1 laboratory by AAS finish, 1 laboratory by ICP-OES finish and 1 laboratory by ICP-MS finish), 50g (1 laboratory by AAS and 2 laboratories by ICP-MS finish) and 200g (4 laboratories by AAS, 1 laboratory by ICP-OES/AAS finish and 1 laboratory by ICP-MS finish);
- Full ICP-OES and ICP-MS elemental suites by 4-acid ($\text{HNO}_3\text{-HF-HClO}_4\text{-HCl}$) digestion (up to 18 laboratories depending on the element);
- Full ICP-OES and ICP-MS elemental suites by aqua regia digestion (up to 19 laboratories depending on the element).

Instrumental neutron activation analysis for Au on 20 x 85mg subsamples was also undertaken at ANSTO, Lucas Heights to confirm homogeneity (see Table 6 below).

Table 3 shows indicative values including major and trace element characterisation based on two samples analysed at Bureau Veritas in Perth, Western Australia which includes:

- Major oxides by lithium borate fusion with X-ray fluorescence;
- LOI at 1000°C by thermogravimetric analyser;
- Total Carbon and Sulphur by infrared combustion furnace;
- Trace elements by laser ablation (on the fused bead) with ICP-MS finish.

For the round robin program twenty 3kg test units were taken at predetermined intervals during the bagging stage, immediately following homogenisation and are considered representative of the entire prepared batch. Six pulp samples were submitted to each laboratory for analysis (the weight provided depended on whether the laboratory was anticipated to undertake assays by gold cyanide leach). The samples received by each laboratory were obtained by taking a sample from six different 3kg test units to maximise representation. The 20 individual INAA results upon which much of the homogeneity evaluation is based, included paired 10g samples taken from 10 different sampling units. This format enabled a nested ANOVA treatment of the INAA results to evaluate homogeneity (see 'Homogeneity Evaluation' section below).

STATISTICAL ANALYSIS

Certified Values and their uncertainty intervals (Tables 1 and 2) have been determined for each analyte following removal of individual, laboratory dataset (batch) and 3SD outliers (single iteration).

For individual outliers within a laboratory batch the z-score test is used in combination with a second method that determines the per cent deviation of the individual value from the

batch median. Outliers in general are selected on the basis of z-scores > 2.5 and with per cent deviations (i) > 3 and (ii) more than three times the average absolute per cent deviation for the batch. Each laboratory data set mean is tested for outlying status based on z-score discrimination and rejected if > 2.5 . After individual and laboratory data set (batch) outliers have been eliminated a non-iterative 3 standard deviation filter is applied, with those values lying outside this window also relegated to outlying status. However, while statistics are taken into account, the exercise of a statistician's prerogative plays a significant role in identifying outliers.

Certified Values are the means of accepted laboratory means after outlier filtering and are the present best estimate of the true value. The INAA data (see Table 6) is omitted from determination of the certified value for Au and is used solely for the calculation of Tolerance Limits and homogeneity evaluation.

95% Expanded Uncertainty provides a 95% probability that the true value of the analyte under consideration lies between the upper and lower limits and is calculated according to the method outlined in ISO 98-3:2008 [6]. All known or suspected sources of bias have been investigated or taken into account.

Indicative (uncertified) values (Table 3) are present where the number of laboratories reporting a particular analyte is insufficient (< 5) to support certification or where interlaboratory consensus is poor.

Homogeneity Evaluation

The tolerance limits (ISO 16269:2014) shown in Tables 1 and 2 were determined using an analysis of precision errors method and are considered a conservative estimate of true homogeneity. The meaning of tolerance limits may be illustrated for copper by 4-acid digestion, where 99% of the time ($1-\alpha=0.99$) at least 95% of subsamples ($p=0.95$) will have concentrations lying between 23.2 and 25.1 ppm. Put more precisely, this means that if the same number of subsamples were taken and analysed in the same manner repeatedly, 99% of the tolerance intervals so constructed would cover at least 95% of the total population, and 1% of the tolerance intervals would cover less than 95% of the total population (ISO Guide 35). ***Please note that tolerance limits pertain to the homogeneity of the CRM only and should not be used as control limits for laboratory performance.***

Table 6 below shows the gold INAA data determined on 20 x 85mg subsamples of OREAS 250c. An equivalent scaled version of the results is also provided to demonstrate an appreciation of what this data means if 30g fire assays were undertaken without the normal measurement error associated with this methodology. In this instance, the 1RSD of 0.233% calculated for a 30g fire assay sample (4.34% at 85mg weights) confirms the high level of gold homogeneity in OREAS 250c.

The homogeneity of gold has been determined by INAA at ANSTO using the reduced analytical subsample method which utilises the known relationship between standard deviation and analytical subsample weight (Ingamells and Switzer, 1973 [2]). In this approach the sample aliquot is substantially reduced to a point where most of the variability in replicate assays should be due to inhomogeneity of the reference material and measurement error becomes negligible.

Table 6. Neutron Activation Analysis of Au (in ppm) on 20 x 85mg subsamples and showing the equivalent results scaled to a 30g sample mass typical of fire assay determination.

Replicate No	Au 85mg actual	Au 30g equivalent*
1	0.315	0.329
2	0.316	0.329
3	0.338	0.330
4	0.351	0.331
5	0.318	0.329
6	0.310	0.328
7	0.366	0.331
8	0.324	0.329
9	0.327	0.329
10	0.330	0.329
11	0.317	0.329
12	0.330	0.329
13	0.322	0.329
14	0.330	0.329
15	0.330	0.329
16	0.323	0.329
17	0.333	0.330
18	0.317	0.329
19	0.337	0.330
20	0.35	0.331
Mean	0.329	0.329
Median	0.328	0.329
Std Dev.	0.014	0.001
Rel.Std.Dev.	4.34%	0.233%

*Results calculated for a 30g equivalent sample mass using the formula: $x^{30g Eq} = \frac{(x^{INAA} - \bar{X}) \times RSD@30g}{RSD@85mg} + \bar{X}$
 where $x^{30g Eq}$ = equivalent result calculated for a 30g sample mass
 (x^{INAA}) = raw INAA result at 85mg
 \bar{X} = mean of 85mg INAA results

The homogeneity of OREAS 250c has also been evaluated in an Analysis of Variance (ANOVA) of the INAA data. The 20 samples were comprised of paired samples from each of 10 sampling lot intervals (representative of the prepared batch) and were randomised prior to assigning sample numbers. The duplicate samples enabled an ANOVA by comparison of within- and between-unit variances across the 10 pairs. The purpose of the ANOVA is to test that no statistically significant difference exists in the variance between units to that of the variance within units. This allows an assessment of homogeneity across the entire prepared batch of OREAS 250c. The test was performed using the following parameters:

- Gold INAA – 20 results (1 laboratory providing duplicate analyses on 10 samples where each sample can be viewed as a ‘unit’);
- Null Hypothesis, H_0 : Between-unit variance is no greater than within-unit variance (reject H_0 if p -value < 0.05);
- Alternative Hypothesis, H_1 : Between-unit variance is greater than within-unit variance.

The data was not filtered for outliers prior to the calculation of the p -value. This process derived a p -value of 0.038, a statistically significant result so the Null Hypothesis is rejected. Any statistically significant result must be checked to determine whether it is also technically significant. A study of the ANOVA components shows that the variation between units is well constrained (between-units RSD at 30g is 0.21%) meaning the magnitude of the effect of between-unit variation is negligible compared to typical (~5%) measurement error associated with the fire assay method.

This isolated case is most likely a false positive (a significant difference is detected where, in reality, none exists) that spuriously manifested from two of the sampling lot pairs falling very close together by chance. The INAA data received also included equally valid ANOVA results for As and La whereby no such significant p -values were detected ($p = 0.449$ and 0.958 , respectively). There is no other supporting evidence to suspect greater between-unit variance compared with within-unit variance. The null hypothesis is therefore retained.

It is important to note that ANOVA is not an absolute measure of homogeneity. Rather, it establishes whether or not the analytes are distributed in a similar manner throughout the packaging run of OREAS 250c and whether the variance between two subsamples from the same unit is statistically distinguishable from the variance of two subsamples taken from any two separate units. A reference material therefore can possess poor absolute homogeneity yet still pass a relative homogeneity (ANOVA) test if the within-unit heterogeneity is large and similar across all units.

Based on the statistical analysis of ANOVA and the results of the interlaboratory certification program, it can be concluded that OREAS 250c is fit-for-purpose as a certified reference material (see 'Intended Use' section below).

PERFORMANCE GATES

The standard deviations (SD's) intervals reported in Table 7 provide an indication of a level of performance that might reasonably be expected from a laboratory being monitored by this CRM in a QA/QC program. They take into account errors attributable to measurement uncertainty and CRM variability. For an effective CRM the contribution of the latter should be negligible in comparison to measurement errors. The Standard Deviation values include all sources of measurement uncertainty: between-lab variance, within-run variance (precision errors) and CRM variability.

In the application of SD's in monitoring performance it is important to note that not all laboratories function at the same level of proficiency and that different methods in use at a particular laboratory have differing levels of precision. Each laboratory has its own inherent SD (for a specific concentration level and analyte-method pair) based on the analytical process and this SD is not directly related to the round robin program (see 'Intended Use' section for more detail). The SD for each analyte's certified value is calculated from the same filtered data set used to determine the certified value, i.e., after removal of all individual, lab dataset (batch) and 3SD outliers (single iteration). These outliers can only be removed after the absolute homogeneity of the CRM has been independently established, i.e., the outliers must be confidently deemed to be analytical rather than arising from inhomogeneity of the CRM.

The standard deviation is then calculated for each analyte from the pooled accepted analyses generated from the certification program.

Table 7 below shows intervals calculated for two and three standard deviations. As a guide these intervals may be regarded as warning or rejection for multiple 2SD outliers, or rejection for individual 3SD outliers in QC monitoring, although their precise application should be at the discretion of the QC manager concerned (also see 'Intended Use' section below). Westgard Rules extend the basics of single-rule QC monitoring using multi-rules (for more information visit www.westgard.com/mltirule.htm). A second method utilises a 5% window calculated directly from the certified value.

Standard deviation is also shown in relative percent for one, two and three relative standard deviations (1RSD, 2RSD and 3RSD) to facilitate an appreciation of the magnitude of these numbers and a comparison with the 5% window. Caution should be exercised when concentration levels approach lower limits of detection of the analytical methods employed as performance gates calculated from standard deviations tend to be excessively wide whereas those determined by the 5% method are too narrow. One approach used at commercial laboratories is to set the acceptance criteria at twice the detection level (DL) \pm 10%.

i.e., Certified Value \pm 10% \pm 2DL [1].

Table 7. Performance Gates for OREAS 250c.

Constituent	Certified Value	Absolute Standard Deviations					Relative Standard Deviations			5% window	
		1SD	2SD Low	2SD High	3SD Low	3SD High	1RSD	2RSD	3RSD	Low	High
Pb Fire Assay											
Au, ppm	0.313	0.015	0.282	0.343	0.266	0.359	4.95%	9.90%	14.84%	0.297	0.328
Aqua Regia Digestion (sample weights 10-50g)											
Au, ppm	0.313	0.018	0.277	0.349	0.259	0.367	5.79%	11.59%	17.38%	0.297	0.329
Cyanide Leach											
Au, ppm	0.295	0.016	0.263	0.328	0.246	0.344	5.51%	11.02%	16.53%	0.280	0.310
PhotonAssay											
Au, ppm	0.311	0.027	0.256	0.366	0.229	0.393	8.79%	17.58%	26.36%	0.295	0.326
4-Acid Digestion											
Ag, ppm	< 0.1	IND	IND	IND	IND	IND	IND	IND	IND	IND	IND
Al, wt. %	5.18	0.101	4.97	5.38	4.87	5.48	1.95%	3.90%	5.85%	4.92	5.43
As, ppm	54	2.2	50	59	47	61	4.15%	8.30%	12.46%	51	57
Ba, ppm	354	7	341	368	334	375	1.93%	3.86%	5.79%	337	372
Be, ppm	1.97	0.087	1.80	2.15	1.71	2.24	4.42%	8.85%	13.27%	1.87	2.07
Bi, ppm	0.29	0.03	0.23	0.35	0.20	0.38	10.16%	20.32%	30.47%	0.28	0.31
Ca, wt. %	1.72	0.060	1.60	1.84	1.54	1.90	3.49%	6.98%	10.48%	1.63	1.80
Cd, ppm	0.075	0.018	0.039	0.110	0.022	0.127	23.62%	47.25%	70.87%	0.071	0.078
Ce, ppm	70	3.1	64	77	61	80	4.33%	8.67%	13.00%	67	74
Co, ppm	20.6	1.32	17.9	23.2	16.6	24.5	6.43%	12.86%	19.29%	19.5	21.6
Cr, ppm	100	6	88	112	82	119	6.10%	12.20%	18.30%	95	105
Cs, ppm	2.62	0.111	2.39	2.84	2.28	2.95	4.25%	8.50%	12.75%	2.49	2.75
Cu, ppm	24.1	2.18	19.8	28.5	17.6	30.7	9.04%	18.09%	27.13%	22.9	25.3
Dy, ppm	3.67	0.204	3.26	4.07	3.05	4.28	5.57%	11.14%	16.71%	3.48	3.85
Er, ppm	1.59	0.107	1.38	1.81	1.27	1.92	6.72%	13.44%	20.16%	1.51	1.67
Eu, ppm	1.60	0.060	1.48	1.72	1.42	1.78	3.74%	7.49%	11.23%	1.52	1.68
Fe, wt. %	4.04	0.136	3.76	4.31	3.63	4.45	3.38%	6.76%	10.14%	3.83	4.24
Ga, ppm	15.3	1.12	13.1	17.6	12.0	18.7	7.33%	14.66%	21.99%	14.6	16.1

SI unit equivalents: ppm (parts per million; 1×10^6) \equiv mg/kg; wt. % (weight per cent) \equiv % (mass fraction).

Note 1: intervals may appear asymmetric due to rounding; IND = indeterminate. Note 2: the number of decimal places quoted does not imply accuracy of the certified value to this level but are given to minimise rounding errors when calculating 2SD and 3SD windows.

Table 7 continued.

Constituent	Certified Value	Absolute Standard Deviations					Relative Standard Deviations			5% window	
		1SD	2SD Low	2SD High	3SD Low	3SD High	1RSD	2RSD	3RSD	Low	High
4-Acid Digestion continued											
Gd, ppm	5.20	0.256	4.69	5.72	4.44	5.97	4.93%	9.86%	14.78%	4.94	5.46
Hf, ppm	4.34	0.289	3.77	4.92	3.48	5.21	6.64%	13.29%	19.93%	4.13	4.56
Ho, ppm	0.63	0.035	0.56	0.70	0.53	0.73	5.51%	11.03%	16.54%	0.60	0.66
In, ppm	0.059	0.004	0.050	0.067	0.045	0.072	7.55%	15.10%	22.64%	0.056	0.061
K, wt. %	1.32	0.052	1.21	1.42	1.16	1.47	3.93%	7.87%	11.80%	1.25	1.38
La, ppm	38.0	1.74	34.5	41.5	32.8	43.2	4.59%	9.19%	13.78%	36.1	39.9
Li, ppm	21.7	1.07	19.6	23.9	18.5	25.0	4.94%	9.88%	14.83%	20.6	22.8
Lu, ppm	0.18	0.02	0.13	0.22	0.11	0.25	12.86%	25.72%	38.59%	0.17	0.19
Mg, wt. %	1.55	0.050	1.45	1.65	1.41	1.70	3.20%	6.40%	9.60%	1.48	1.63
Mn, wt. %	0.049	0.001	0.046	0.052	0.045	0.054	3.02%	6.05%	9.07%	0.047	0.052
Mo, ppm	1.96	0.167	1.63	2.30	1.46	2.47	8.50%	17.00%	25.49%	1.87	2.06
Na, wt. %	1.07	0.023	1.03	1.12	1.00	1.14	2.18%	4.36%	6.54%	1.02	1.13
Nb, ppm	31.8	1.77	28.2	35.3	26.5	37.1	5.57%	11.15%	16.72%	30.2	33.4
Nd, ppm	32.1	1.70	28.7	35.5	27.0	37.2	5.30%	10.60%	15.90%	30.5	33.7
Ni, ppm	83	4.0	75	91	71	95	4.83%	9.66%	14.49%	79	87
P, wt. %	0.095	0.003	0.088	0.102	0.085	0.106	3.57%	7.14%	10.71%	0.091	0.100
Pb, ppm	9.92	0.915	8.09	11.75	7.18	12.66	9.22%	18.44%	27.65%	9.42	10.42
Pr, ppm	8.51	0.465	7.58	9.44	7.12	9.91	5.46%	10.92%	16.38%	8.09	8.94
Rb, ppm	65	3.8	57	72	53	76	5.87%	11.74%	17.61%	61	68
Re, ppm	< 0.002	IND	IND	IND	IND	IND	IND	IND	IND	IND	IND
S, wt. %	0.027	0.004	0.020	0.035	0.016	0.038	13.71%	27.41%	41.12%	0.026	0.029
Sb, ppm	0.88	0.068	0.75	1.02	0.68	1.09	7.70%	15.40%	23.11%	0.84	0.93
Sc, ppm	9.89	0.624	8.64	11.14	8.02	11.76	6.31%	12.61%	18.92%	9.40	10.39
Sm, ppm	6.27	0.336	5.60	6.94	5.26	7.28	5.36%	10.72%	16.07%	5.96	6.59
Sn, ppm	3.59	0.166	3.26	3.92	3.09	4.09	4.61%	9.23%	13.84%	3.41	3.77
Sr, ppm	311	8	296	327	288	335	2.48%	4.96%	7.45%	296	327
Ta, ppm	2.07	0.175	1.73	2.42	1.55	2.60	8.43%	16.86%	25.28%	1.97	2.18
Tb, ppm	0.69	0.018	0.65	0.72	0.63	0.74	2.64%	5.28%	7.92%	0.65	0.72
Th, ppm	8.68	0.344	7.99	9.36	7.64	9.71	3.96%	7.92%	11.88%	8.24	9.11
Ti, wt. %	0.567	0.014	0.539	0.595	0.525	0.609	2.47%	4.94%	7.41%	0.539	0.595
Tl, ppm	0.35	0.018	0.32	0.39	0.30	0.41	5.00%	10.00%	15.00%	0.34	0.37
Tm, ppm	0.20	0.013	0.18	0.23	0.16	0.24	6.26%	12.51%	18.77%	0.19	0.21
U, ppm	1.62	0.107	1.41	1.84	1.30	1.94	6.59%	13.18%	19.77%	1.54	1.70
V, ppm	70	2.6	65	75	62	78	3.69%	7.38%	11.07%	66	73
W, ppm	2.05	0.137	1.78	2.33	1.64	2.46	6.66%	13.32%	19.97%	1.95	2.16
Y, ppm	15.7	0.96	13.7	17.6	12.8	18.6	6.14%	12.29%	18.43%	14.9	16.4
Yb, ppm	1.28	0.063	1.16	1.41	1.09	1.47	4.90%	9.79%	14.69%	1.22	1.35
Zn, ppm	75	3.3	68	82	65	85	4.42%	8.84%	13.26%	71	79
Zr, ppm	188	10	167	208	157	219	5.47%	10.94%	16.41%	178	197
Aqua Regia Digestion											
Ag, ppm	0.079	0.008	0.063	0.094	0.056	0.102	9.73%	19.45%	29.18%	0.075	0.083
Al, wt. %	1.64	0.122	1.40	1.89	1.28	2.01	7.41%	14.83%	22.24%	1.56	1.72
As, ppm	43.4	1.33	40.8	46.1	39.4	47.4	3.07%	6.14%	9.21%	41.3	45.6
B, ppm	< 10	IND	IND	IND	IND	IND	IND	IND	IND	IND	IND

SI unit equivalents: ppm (parts per million; 1×10^{-6}) \equiv mg/kg; wt. % (weight per cent) \equiv % (mass fraction).

Note 1: intervals may appear asymmetric due to rounding; IND = indeterminate. Note 2: the number of decimal places quoted does not imply accuracy of the certified value to this level but are given to minimise rounding errors when calculating 2SD and 3SD windows.

Table 7 continued.

Constituent	Certified Value	Absolute Standard Deviations					Relative Standard Deviations			5% window	
		1SD	2SD Low	2SD High	3SD Low	3SD High	1RSD	2RSD	3RSD	Low	High
Aqua Regia Digestion continued											
Ba, ppm	95	5.0	85	104	80	109	5.25%	10.51%	15.76%	90	99
Be, ppm	1.08	0.069	0.94	1.22	0.87	1.29	6.43%	12.86%	19.30%	1.02	1.13
Bi, ppm	0.26	0.021	0.22	0.31	0.20	0.33	7.85%	15.70%	23.55%	0.25	0.28
Ca, wt. %	0.524	0.041	0.442	0.606	0.401	0.647	7.81%	15.63%	23.44%	0.498	0.551
Cd, ppm	0.067	0.023	0.021	0.114	0.000	0.137	34.70%	69.39%	104.09%	0.064	0.071
Ce, ppm	40.9	1.87	37.1	44.6	35.3	46.5	4.59%	9.17%	13.76%	38.8	42.9
Co, ppm	15.9	0.71	14.5	17.4	13.8	18.1	4.49%	8.97%	13.46%	15.1	16.7
Cr, ppm	45.9	5.4	35.0	56.8	29.6	62.2	11.83%	23.66%	35.49%	43.6	48.2
Cs, ppm	1.13	0.093	0.94	1.31	0.85	1.41	8.22%	16.44%	24.66%	1.07	1.18
Cu, ppm	18.9	2.2	14.5	23.2	12.4	25.4	11.51%	23.02%	34.53%	17.9	19.8
Fe, wt. %	3.01	0.263	2.48	3.53	2.22	3.80	8.74%	17.48%	26.21%	2.86	3.16
Ga, ppm	5.97	0.521	4.93	7.01	4.41	7.53	8.73%	17.45%	26.18%	5.67	6.27
Ge, ppm	0.100	0.017	0.066	0.133	0.049	0.150	16.97%	33.94%	50.91%	0.095	0.105
In, ppm	0.025	0.003	0.019	0.032	0.015	0.036	13.46%	26.91%	40.37%	0.024	0.027
K, wt. %	0.335	0.025	0.285	0.385	0.260	0.409	7.41%	14.82%	22.24%	0.318	0.352
La, ppm	23.2	0.97	21.2	25.1	20.3	26.1	4.17%	8.34%	12.50%	22.0	24.3
Li, ppm	7.61	0.628	6.35	8.86	5.72	9.49	8.25%	16.50%	24.75%	7.23	7.99
Mg, wt. %	0.989	0.069	0.851	1.126	0.782	1.195	6.95%	13.90%	20.85%	0.939	1.038
Mn, wt. %	0.033	0.002	0.029	0.038	0.026	0.040	7.00%	14.01%	21.01%	0.032	0.035
Mo, ppm	1.32	0.118	1.09	1.56	0.97	1.68	8.95%	17.90%	26.85%	1.26	1.39
Na, wt. %	0.351	0.022	0.307	0.394	0.285	0.416	6.26%	12.53%	18.79%	0.333	0.368
Nb, ppm	0.59	0.09	0.40	0.77	0.30	0.87	16.03%	32.06%	48.10%	0.56	0.62
Ni, ppm	73	3.4	66	79	62	83	4.73%	9.45%	14.18%	69	76
P, wt. %	0.058	0.004	0.051	0.066	0.047	0.070	6.50%	13.00%	19.50%	0.056	0.061
Pb, ppm	7.08	0.556	5.97	8.19	5.41	8.75	7.85%	15.69%	23.54%	6.73	7.44
Rb, ppm	19.1	1.16	16.8	21.4	15.6	22.5	6.07%	12.14%	18.20%	18.1	20.0
Re, ppm	< 0.001	IND	IND	IND	IND	IND	IND	IND	IND	IND	IND
S, wt. %	0.024	0.004	0.016	0.033	0.012	0.037	17.46%	34.93%	52.39%	0.023	0.026
Sb, ppm	0.45	0.06	0.33	0.56	0.27	0.62	13.05%	26.11%	39.16%	0.42	0.47
Sc, ppm	3.74	0.314	3.11	4.37	2.79	4.68	8.41%	16.82%	25.23%	3.55	3.92
Sn, ppm	1.24	0.13	0.98	1.50	0.84	1.63	10.58%	21.15%	31.73%	1.17	1.30
Sr, ppm	63	2.6	58	69	56	71	4.08%	8.17%	12.25%	60	67
Ta, ppm	< 0.01	IND	IND	IND	IND	IND	IND	IND	IND	IND	IND
Tb, ppm	0.45	0.034	0.38	0.52	0.35	0.55	7.62%	15.25%	22.87%	0.43	0.47
Th, ppm	5.85	0.247	5.36	6.34	5.11	6.59	4.21%	8.43%	12.64%	5.56	6.14
Ti, wt. %	0.162	0.013	0.137	0.188	0.124	0.201	7.88%	15.75%	23.63%	0.154	0.170
Tl, ppm	0.14	0.02	0.10	0.17	0.09	0.18	11.67%	23.35%	35.02%	0.13	0.14
U, ppm	0.97	0.079	0.81	1.12	0.73	1.20	8.15%	16.31%	24.46%	0.92	1.01
V, ppm	33.7	1.83	30.0	37.4	28.2	39.2	5.43%	10.86%	16.29%	32.0	35.4
W, ppm	0.38	0.038	0.31	0.46	0.27	0.50	9.99%	19.98%	29.97%	0.36	0.40
Y, ppm	9.10	0.386	8.33	9.87	7.94	10.26	4.24%	8.48%	12.71%	8.65	9.56
Yb, ppm	0.63	0.040	0.55	0.71	0.51	0.75	6.39%	12.78%	19.16%	0.60	0.66
Zn, ppm	46.4	4.7	37.0	55.8	32.3	60.5	10.13%	20.26%	30.39%	44.1	48.7
Zr, ppm	37.1	8.6	19.8	54.4	11.2	63.0	23.29%	46.57%	69.86%	35.3	39.0

SI unit equivalents: ppm (parts per million; 1×10^{-6}) \equiv mg/kg; wt. % (weight per cent) \equiv % (mass fraction).

Note 1: intervals may appear asymmetric due to rounding; IND = indeterminate.

Note 2: the number of decimal places quoted does not imply accuracy of the certified value to this level but are given to minimise rounding errors when calculating 2SD and 3SD windows.

PARTICIPATING LABORATORIES

1. Actlabs, Ancaster, Ontario, Canada
2. Alex Stewart International, Mendoza, Argentina
3. ALS, Canning Vale, WA, Australia
4. ALS, Johannesburg, South Africa
5. ALS, Kalgoorlie, WA, Australia
6. ALS, Lima, Peru
7. ALS, Loughrea, Galway, Ireland
8. ALS, Malaga, WA, Australia
9. ALS, Vancouver, BC, Canada
10. American Assay Laboratories, Sparks, Nevada, USA
11. ANSTO, Lucas Heights, NSW, Australia
12. Bureau Veritas Geoanalytical, Perth, WA, Australia
13. CERTIMIN, Lima, Peru
14. Gekko Assay Labs, Ballarat, VIC, Australia
15. Inspectorate (BV), Lima, Peru
16. Intertek Genalysis, Perth, WA, Australia
17. Intertek Tarkwa, Tarkwa, Ghana
18. Intertek Testing Services Philippines, Cupang, Muntinlupa, Philippines
19. MSA ENVAL Laboratories, Yamoussoukro, Côte d'Ivoire
20. MSALABS, Prince George, BC, Canada
21. MSALABS, Val-d'Or, Quebec, Canada
22. MSALABS Bulyanhulu Gold Mine, Bubada, Shinyanga, United Republic of Tanzania
23. MSALABS Kibali Gold Mines, Doko, Haut-Uélé, Congo, Democratic Republic of the (Zaire)
24. On Site Laboratory Services, Bendigo, VIC, Australia
25. PT Geoservices Ltd, Cikarang, Jakarta Raya, Indonesia
26. PT Intertek Utama Services, Jakarta Timur, DKI Jakarta, Indonesia
27. Ravenswood Gold, Ravenswood, QLD, Australia
28. Reminex Centre de Recherche, Marrakesh, Marrakesh-Safi, Morocco
29. SGS, Randfontein, Gauteng, South Africa
30. SGS Australia Mineral Services, Perth, WA, Australia
31. SGS Tarkwa, Tarkwa, Western Region, Ghana
32. Shiva Analyticals Ltd, Bangalore North, Karnataka, India
33. Stewart Assay & Environmental Laboratories LLC, Kara-Balta, Chüy, Kyrgyzstan

Please note: To preserve anonymity, the above numbered alphabetical list of participating laboratories does not correspond with the Lab ID numbering on the scatter plots below.

PREPARER AND SUPPLIER

Certified reference material OREAS 250c is prepared, certified and supplied by:



ORE Research & Exploration Pty Ltd
37A Hosie Street
Bayswater North VIC 3153
AUSTRALIA

Tel: +613-9729 0333
Fax: +613-9729 8338
Web: www.oreas.com
Email: info@ore.com.au

Figure 1. Au by Pb Fire Assay in OREAS 250c

SPC.1680.RR-250~c260~b.OREAS 250c.2.Fire Assay.Au.Lab.231010.163432.SN

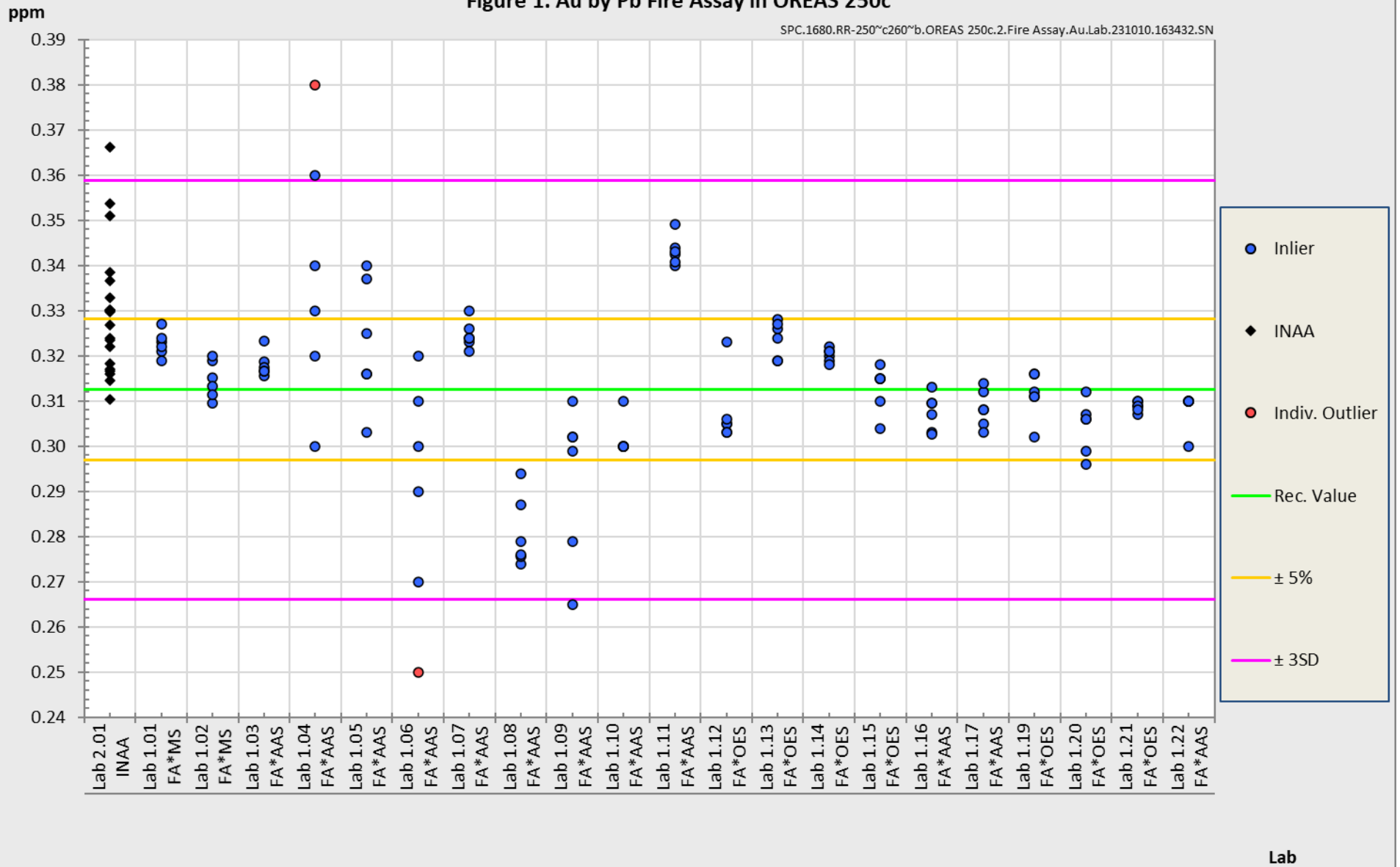


Figure 2. Au by Aqua Regia Digestion (sample weights 10-50g) in OREAS 250c

SPC.1680.RR-250~c260~b.OREAS 250c.2.AR Digest 10-50g.Au.Lab.231010.171413.SN

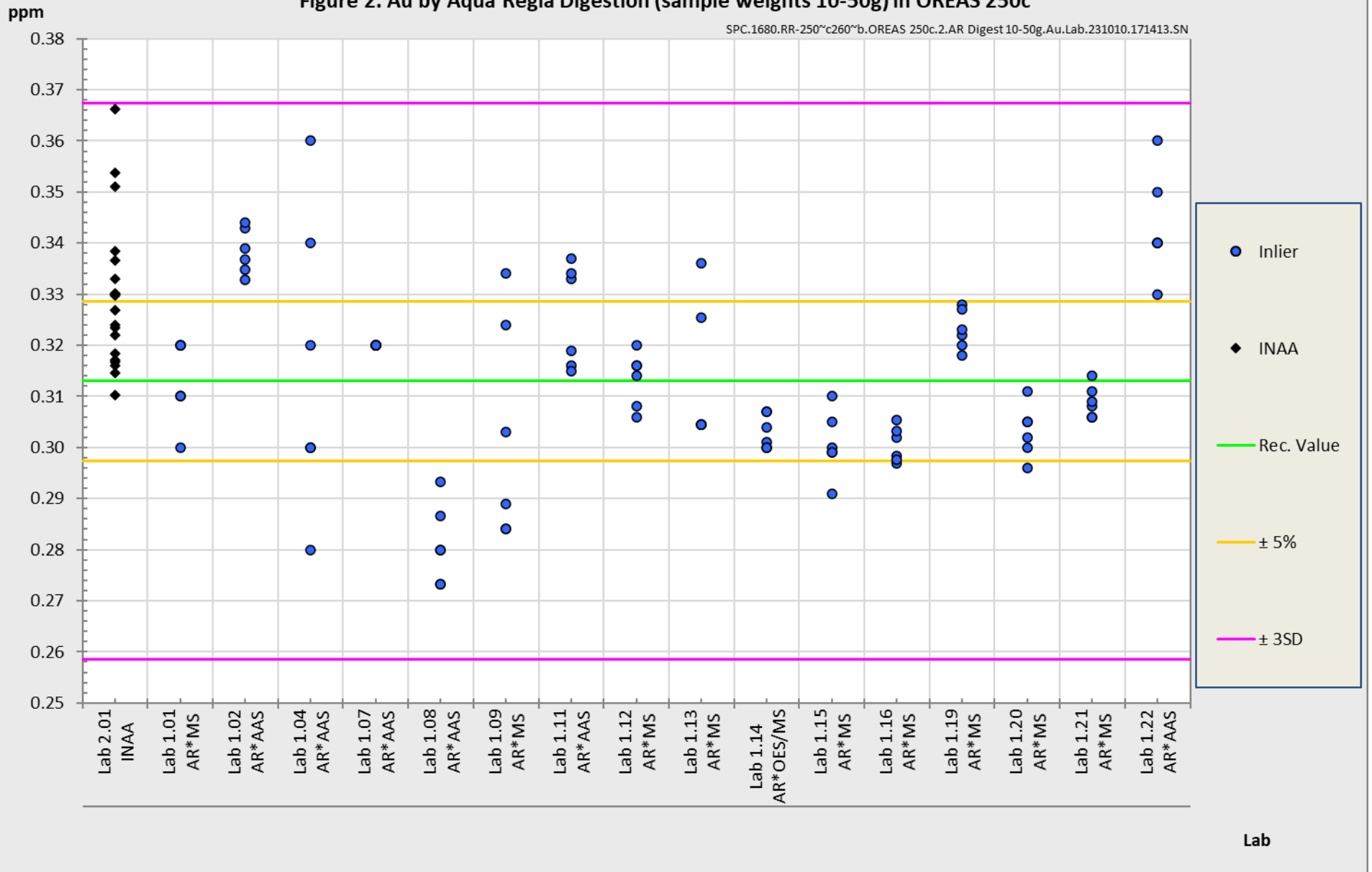


Figure 3. Au by Cyanide Leach in OREAS 250c

SPC.1680.RR-250~c260~b.OREAS 250c.2.CNL.Au.Lab.231010.172220.SN


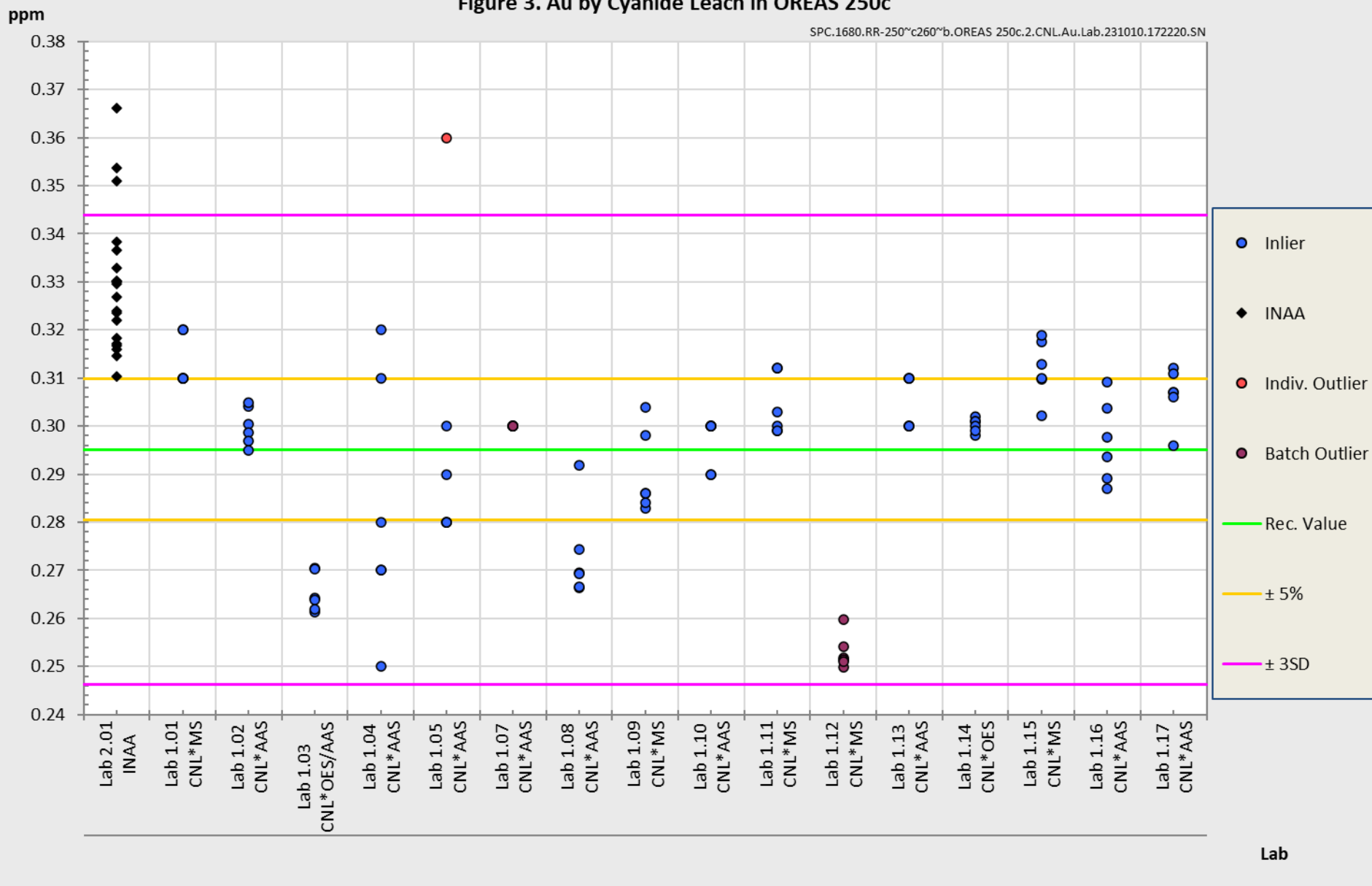


Figure 4. Au by PhotonAssay in OREAS 250c

SPC.1680.RR-250~c260~b.OREAS 250c.3.PA.Au.Lab.240228.130944.SN

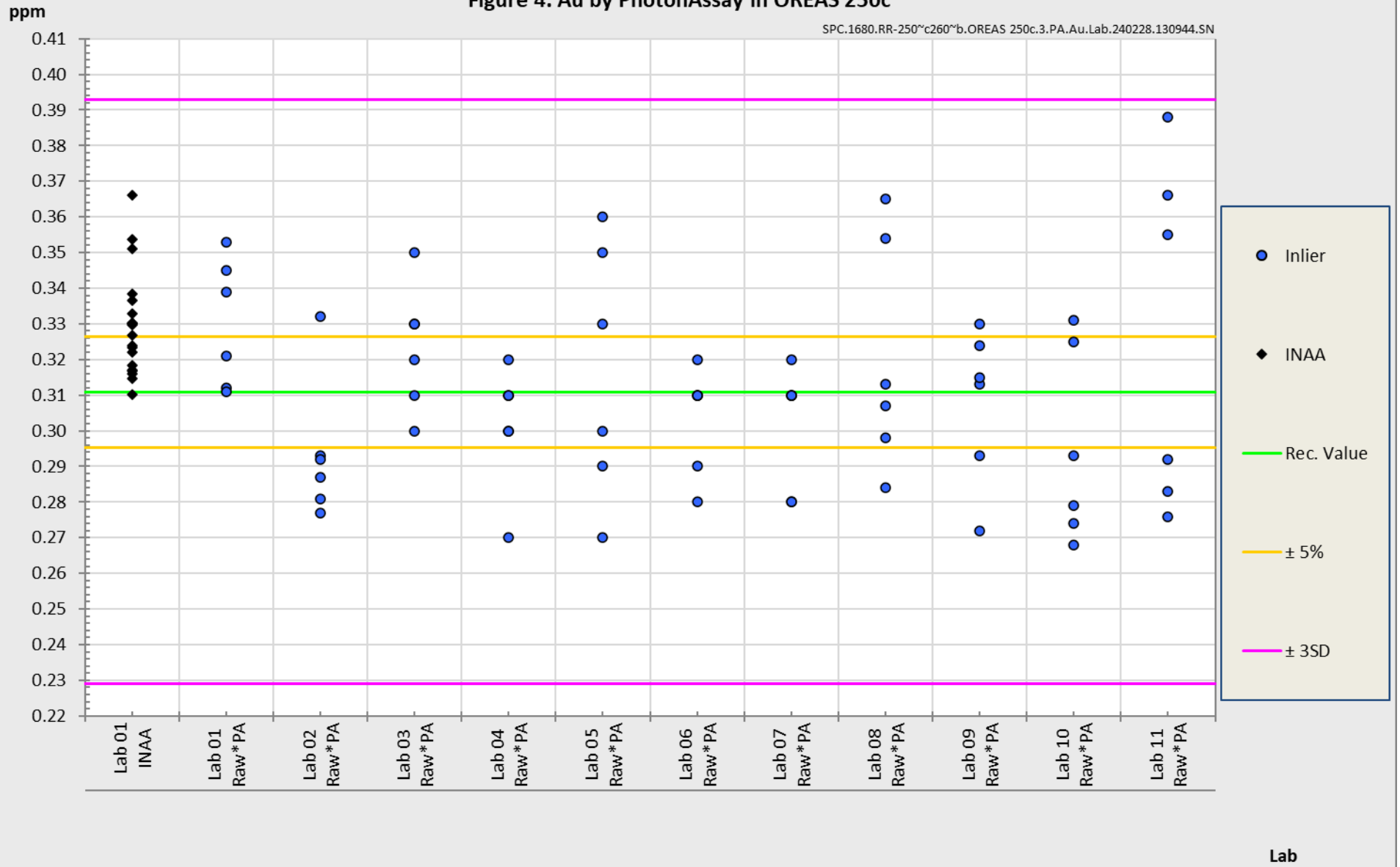

METROLOGICAL TRACEABILITY

The interlaboratory results that underpin the certified values are metrologically traceable to the international measurement scale (SI) of mass (either as a % mass fraction or as milligrams per kilogram (mg/kg)). In line with popular use, all data within tables in this certificate are expressed as the mass fraction in either weight percent (wt.%) or parts per million (ppm).

The analytical samples sent to participating laboratories were selected in a manner to be representative of the entire prepared batch of CRM. This 'representivity' was maintained in each submitted laboratory sample batch and ensures the user that the data is traceable from sample selection through to the analytical results. The systematic sampling method was chosen due to the low risk of overlooking repetitive effects or trends in the batch due to the way the CRM was processed. In line with ISO 17025 [10], each analytical data set received from the participating laboratories has been validated by its assayer through the inclusion of internal reference materials and QC checks during and post analysis.

The participating laboratories were chosen on the basis of their competence (from past performance in interlaboratory programs undertaken by ORE Pty Ltd) for a particular analytical method, analyte or analyte suite and sample matrix. These laboratories are accredited to ISO 17025 for Au by fire assay (Table 1). The other operationally defined measurands characterised in this certificate (Table 2) are derived from data procured mostly from ISO 17025 accredited laboratories. The certified values presented in this report are calculated from the means of accepted data following robust technical and statistical analysis as detailed in this report.

Guide ISO/TR 16476:2016, section 5.3.1 describes metrological traceability in reference materials as it pertains to the transformation of the measurand. In this section it states, *"Although the determination of the property value itself can be made traceable to appropriate units through, for example, calibration of the measurement equipment used, steps like the transformation of the sample from one physical (chemical) state to another cannot. Such transformations may only be compared with a reference (when available), or among themselves. For some transformations, reference methods have been defined and may be used in certification projects to evaluate the uncertainty associated with such a transformation. In other cases, only a comparison among different laboratories using the same procedure is possible. In this case, it is impossible to demonstrate absence of method bias; therefore, the result is an operationally defined measurand (ISO Guide 35:2017, 9.2.4c)."* Certification takes place on the basis of agreement among operationally defined, independent measurement results.

COMMUTABILITY

The measurements of the results that underlie the certified values contained in this report were undertaken by methods involving pre-treatment (fusion/digestion) of the sample. This served to reduce the sample to a simple and well understood form permitting calibration using simple solutions of the CRM. Due to these methods being well understood and highly effective, commutability is not an issue for this CRM. All OREAS CRMs are sourced from natural ore minerals meaning they will display similar behaviour as routine 'field' samples in the relevant measurement process. Care should be taken to ensure 'matrix matching' as close as practically achievable. The matrix and mineralisation style of the CRM is described

in the 'Source Material' section and users should select appropriate CRMs matching these attributes to the field samples being analysed.

INTENDED USE

OREAS 250c is intended to cover all activities needed to produce a measurement result. This includes extraction, possible separation steps and the actual measurement process (the signal producing step). OREAS 250c may be used to calibrate the entire procedure by producing a pure substance CRM transformed into a calibration solution.

OREAS 250c is intended for the following uses:

- For the monitoring of laboratory performance in the analysis of analytes reported in Tables 1 and 2 in geological samples;
- For the verification of analytical methods for analytes reported in Tables 1 and 2;
- For the calibration of instruments used in the determination of the concentration of analytes reported in Tables 1 and 2. When a value provided in this certificate is used to calibrate a measurement process, the uncertainty associated with that value should be appropriately propagated into the user's uncertainty calculation. Users can determine an approximation of the standard uncertainty by calculating one fourth of the width of the Expanded Uncertainty interval given in this certificate (Expanded Uncertainty intervals are provided in Tables 1 and 2).

MINIMUM SAMPLE SIZE

To relate analytical determinations to the values in this certificate, the minimum mass of sample used should match the typical mass that the laboratories used in the interlaboratory (round robin) certification program. This means that different minimum sample masses should be used depending on the operationally defined methodology as follows:

- Au by fire assay: $\geq 15\text{g}$;
- Au by aqua regia digestion: $\geq 10\text{g}$;
- Au by cyanide leach: $\geq 5\text{g}$;
- Au by PhotonAssay: $\sim 350\text{g}$;
- 4-acid digestion with ICP-OES and/or MS finish: $\geq 0.25\text{g}$;
- Aqua regia digestion with ICP-OES and/or MS finish: $\geq 0.5\text{g}$.

PERIOD OF VALIDITY & STORAGE INSTRUCTIONS

The certification of OREAS 250c remains valid, within the specified measurement uncertainties, until July 2038, provided the CRM is handled and stored in accordance with the instructions given below. This certification is nullified if the CRM is any way changed or contaminated.

Store in a clean and cool dry place away from direct sunlight.

Long-term stability will be monitored at appropriate intervals and purchasers notified if any changes are observed. The period of validity may well be indefinite and will be reassessed prior to expiry with the aim of extending the validity if possible.

Single-use sachets

Following analysis, it is the manufacturer's expectation that any remaining material is discarded unless the sachet is promptly resealed. It is the user's responsibility to prevent contamination and minimise exposure to the atmosphere.

Repeat-use packaging (e.g., 500g unit)

After taking a subsample, users should replace the lid of the jar promptly and securely to prevent accidental spills and airborne contamination. OREAS 250c contains a non-hygroscopic* matrix with an indicative value for moisture provided to enable users to check for changes to stored material by determining moisture in the user's laboratory and comparing the result to the value in Table 4 in this certificate.

The stability of the CRM in regard to oxidation from the breakdown of sulphide minerals to sulphates is negligible given its low sulphur concentration (0.01 wt.% S).

*A non-hygroscopic matrix means exposure to atmospheres significantly different, in terms of temperature and humidity, from the climate during manufacturing should have negligible impact on the precision of results. Hygroscopic moisture is the amount of adsorbed moisture (weakly held H₂O- molecules on the surface of exposed material) following exposure to the local atmosphere. Usually, equilibration of material to the local atmosphere will only occur if the material is spread into a thin (~2mm thick) layer and left exposed for a period of 2 hours.

INSTRUCTIONS FOR HANDLING & CORRECT USE

Pre-homogenisation of the CRM prior to subsampling and analysis is not necessary as there is no particle segregation under transport [13].

Fine powders pose a risk to eyes and lungs and therefore standard precautions including the use of safety glasses and dust masks are advised.

QC monitoring using multiples of the Standard Deviation (SD)

In the application of SD's in monitoring performance it is important to note that not all laboratories function at the same level of proficiency and that different methods in use at a particular laboratory have differing levels of precision. Each laboratory has its own inherent SD (for a specific concentration level and analyte-method pair) based on the analytical process and this SD is not directly related to the round robin program.

The majority of data generated in the round robin program was produced by a selection of world class laboratories. The SD's thus generated are more constrained than those that would be produced across a randomly selected group of laboratories. To produce more generally achievable SD's the 'pooled' SD's provided in this report include interlaboratory bias. This 'one size fits all' approach may require revision at the discretion of the QC manager concerned following careful scrutiny of QC control charts.

The performance gates shown in Table 7 are intended only to be used as a preliminary guide as to what a laboratory may be able to achieve. Over a period of time monitoring your own laboratory's data for this CRM, SD's should be calculated directly from your own laboratory's process. This will enable you to establish more specific performance gates that are fit for purpose for your application as well as the ability to monitor bias. If your long-term trend analysis shows an average value that is within the 95% expanded uncertainty then generally there is no cause for concern in regard to bias.

For use with the aqua regia digestion method

It is important to note that in the analytical industry there is no standardisation of the aqua regia digestion process. This method is a partial empirical digest and differences in recoveries for various analytes are commonplace. These are caused by variations in the digest conditions and can include the ratio of nitric to hydrochloric acids, acid strength, temperatures, leach times and secondary digestions. Recoveries for sulphide-hosted base metal sulphides approach total values, however, other analytes, in particular the lithophile elements, show greater sensitivity to method parameters. This can result in lack of consensus in an inter-laboratory certification program for these elements.

The approach applied here is to report certified values in those instances where reasonable agreement exists amongst a majority of participating laboratories. The results of specific laboratories may differ significantly from the certified values, but will, nonetheless, be valid and reproducible in the context of the specifics of the aqua regia method in use. Users of this reference material should, therefore, be mindful of this limitation when applying the certified values in a quality control program.

LEGAL NOTICE

Ore Research & Exploration Pty Ltd has prepared and statistically evaluated the property values of this reference material to the best of its ability. The Purchaser by receipt hereof releases and indemnifies Ore Research & Exploration Pty Ltd from and against all liability and costs arising from the use of this material and information.

© COPYRIGHT Ore Research & Exploration Pty Ltd.
Unauthorised copying, reproduction, storage or dissemination is prohibited.

DOCUMENT HISTORY

Revision No.	Date	Changes applied
2	18 th March, 2024	Added Au by PhotonAssay certification.
1	15 th December, 2023	Added SQ-XRD mineralogical data to Table 5.
0	16 th October, 2023	First publication.

CERTIFYING OFFICER



18th March, 2024

Craig Hamlyn (B.Sc. Hons - Geology), Technical Manager - ORE P/L

QMS CERTIFICATION

ORE Pty Ltd is accredited for compliance with ISO 17034:2016.



ORE Pty Ltd is ISO 9001:2015 certified by Lloyd's Register Quality Assurance Ltd for its quality management system including development, manufacturing, certification and supply of CRMs.



REFERENCES

- [1] Govett, G.J.S. (1983). Handbook of Exploration Geochemistry, Volume 2: Statistics and Data Analysis in Geochemical Prospecting (Variations of accuracy and precision).
- [2] Ingamells, C. O. and Switzer, P. (1973). A Proposed Sampling Constant for Use in Geochemical Analysis, Talanta 20, 547-568.
- [3] ISO Guide 30:2015. Terms and definitions used in connection with reference materials.
- [4] ISO Guide 31:2015. Reference materials – Contents of certificates and labels.
- [5] ISO Guide 35:2017. Certification of reference materials - General and statistical principals.
- [6] ISO Guide 98-3:2008. Guide to the expression of uncertainty in measurement (GUM:1995).
- [7] ISO 16269:2014. Statistical interpretation of data – Part 6: Determination of statistical tolerance intervals.
- [8] ISO/TR 16476:2016, Reference Materials – Establishing and expressing metrological traceability of quantity values assigned to reference materials.
- [9] ISO 17025:2017, General requirements for the competence of testing and calibration laboratories.
- [10] ISO 17034:2016. General requirements for the competence of reference material producers.
- [11] Munsell Rock Color Book (2014). Rock-Color Chart Committee, Geological Society of America (GSA), Minnesota (USA).
- [12] OREAS-BUP-70-09-11: Statistical Analysis - OREAS Evaluation Method.
- [13] OREAS-TN-04-1498: Stability under transport; an experimental study of OREAS CRMs.
- [14] OREAS-TN-05-1674: Long-term storage stability; an experimental study of OREAS CRMs.
- [15] Thompson, A.; Taylor, B.N.; Guide for the Use of the International System of Units (SI); NIST Special Publication 811; U.S. Government Printing Office: Washington, DC (2008); available at: <https://physics.nist.gov/cuu/pdf/sp811.pdf> (accessed Nov 2021).
- [16] Van der Veen AMH and Pauwels, J. (2001). Uncertainty calculations in the certification of reference materials, Accred Qual Assur 6: 290-294.